

V.A.C. Via™ Negative Pressure Wound Therapy System

QUICK REFERENCE GUIDE



V.A.C. Via™ Unterdruckwundtherapie-System Kurzanleitung
 Naslaggids voor V.A.C. Via™ Negative Pressure Wound Therapy System
 Quick Reference Guide du système de traitement des plaies par pression négative V.A.C. Via™
 Guida di riferimento rapido del sistema di terapia a pressione negativa V.A.C. Via™
 Guía de Referencia Rápida del Sistema de Terapia de Presión Negativa para el Tratamiento de Heridas V.A.C. Via™
 Lynvejledning til V.A.C. Via™ Undertryksbehandlingssystem til sår
 Snabbreferensguide för V.A.C. Via™ sårbehandlingssystem med negativt tryck
 Guia de Referência Rápida do Sistema de Terapia de Lesões com Pressão Negativa V.A.C. Via™
 V.A.C. Via™ Negatif Basınc Yara Terapi Sistemi Hızlı Ba vurur Kılavuzu
 Οδηγός γρήγορης αναφοράς για το Σύστημα θεραπείας επούλωσης τραυμάτων με εφαρμογή αρνητικής πίεσης V.A.C. Via™
 V.A.C. Via™-alipaineimuhoitojärjestelmän pikaopas
 Hurtigveiledning for V.A.C. Via™-sårbehandlingssystem med negativt trykk

7-DAY V.A.C.® THERAPY SYSTEM

7-Tage V.A.C.® Therapy-System
 7-daags V.A.C.® Therapy System
 Système V.A.C.® Therapy jours
 Sistema V.A.C.® Therapy settimanale
 Sistema de Terapia V.A.C.® en 7 Días
 7-dages V.A.C.® Terapisystem

7-dagars V.A.C.® Therapy-system
 Sistema de Terapia de 7 dias V.A.C.®
 7-Günlük V.A.C.® Tedavi Sistemi
 Σύστημα Θεραπείας 7 ημερών V.A.C.®
 7 päivän V.A.C.® -hoitojärjestelmä
 7-dagers V.A.C.® -behandlingssystem



V.A.C.Via™ Negative Pressure Wound Therapy System - Quick Reference Guide

The V.A.C.Via™ Therapy Unit is a single use, disposable V.A.C.® Therapy device designed for moderate to low-severity wounds. The V.A.C.Via™ Therapy Unit has a 7-day lifespan and a rechargeable battery. This unit provides negative pressure at either 75mmHg or 125mmHg, and offers selection of Continuous therapy or Dynamic Pressure Control™ therapy. The V.A.C.Via™ Starter Kits include one therapy unit, one carrying case, one AC power cord, one 250cc canister and one Medium Spiral Granufoam Dressing. Dressing kits and canisters designed specifically for the V.A.C.Via™ are provided separately.

CAUTION: Once therapy is on for one continuous hour, the 7-day lifespan begins (one light for each day of therapy remaining) and continues even if unit is turned off.

IMPORTANT: As with any prescription medical device, failure to consult a physician and carefully read and follow all therapy unit and dressing instructions and safety information prior to each use may lead to improper product performance and the potential for serious or fatal injury. Do not adjust therapy unit settings or perform therapy application without directions from/or supervision by the treating physician or caregiver.

This guide is not intended as a comprehensive resource. For additional information concerning proper use of the V.A.C.Via™ Therapy System and important information including Safety Information, Cautions and Warnings, consult the Instructions for Use located in the device packaging, and discuss with your doctor.

Upon receipt, of the V.A.C.Via™ Therapy Unit battery may not be fully charged.


Dressing Application



NOTE: Ensure clinician records the total number of pieces of foam used in the wound and documents on the supplied Foam Quantity Label (attached to the Interface Pad tubing) and in the patient's chart.

NOTE: The use of 3M™ Cavilon™ No Sting Barrier Film is suggested and may be provided with V.A.C.Via™ dressings in some but not all countries, to protect and prepare periwound skin.

Canister Installation/Removal

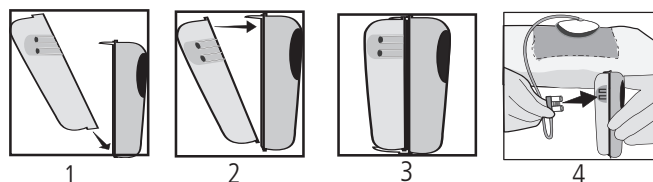


To remove canister, disconnect tubing and depress locking tab on top of canister.

All disposable components of the V.A.C.® Therapy System are for single use only. Re-use of disposable components may result in wound contamination, infection and / or failure of the wound to heal.


Turn Therapy On

Press and Hold center of On Button for Three Seconds (125mmHg LED Will Flash Until Pressure is Reached)

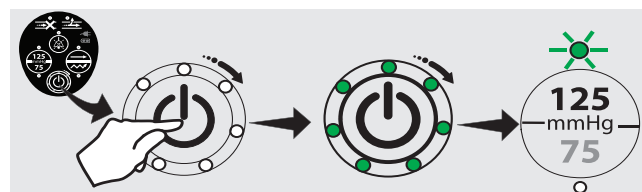

Therapy Settings


Select Pressure Setting 125mmHg (Default) or 75mmHg (Press and Hold center for Three Seconds)

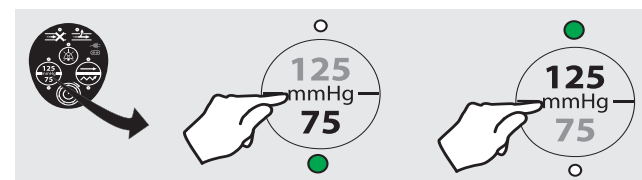

Select Continuous (Default) or Dynamic Pressure Control™ (Press and Hold center for Three Seconds)




Therapy Unit Carrying Case

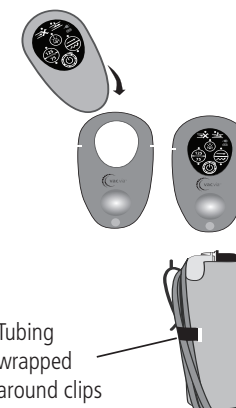

The V.A.C.Via™ Carrying Case has an integrated belt clip and a separate adjustable strap to allow for versatile positioning.

CAUTION: Do not wear strap around neck.

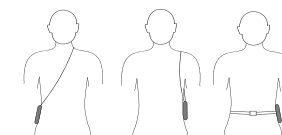

Tubing management - the belt clip and additional clips on each side and at the bottom of the carrying case provide a place where excess tubing may be wrapped and stored.

Instructions for Patient

For important **Safety Information, Warnings, Precautions** and detailed operating Instructions see V.A.C.Via™ **Instructions for Use** (located in the device packaging). Always consult a physician and carefully review device instructions and Safety Information with your clinical care provider.

DAILY USE

The V.A.C.Via™ Therapy System is portable and small enough that it may be worn beneath clothing during normal patient activities as approved by the treating physician.

Sleeping:

- Position therapy unit so that tubing will not become kinked or pinched.
- Ensure therapy unit will not be pulled off a table or fall to the floor during sleep.
- It is recommended to keep therapy unit plugged in and charging while sleeping.

Showering and Bathing:

- Do not use the V.A.C.Via™ Therapy Unit where it can fall or be pulled into a tub, shower or sink.
- Do not reach for a product that has fallen into water. Unplug the unit immediately if plugged into an electrical source. Disconnect unit from dressing and contact treating Physician or Caregiver.
- To shower, clamp tubing and disconnect from therapy unit.
- The clear drape is waterproof; patient may wash or shower with dressing in place and the tube clamped and disconnected from the therapy unit.
- When towel drying avoid disturbing or damaging the dressing.

Cleaning:

- The V.A.C.Via™ Therapy Unit and Carrying Case can be wiped with a damp cloth using a mild household cleaner.

KEEP GUIDE WITH DEVICE

Therapy Life Indicators

Days Remaining of the 7-day Therapy Cycle



CAUTION: Once therapy is on for one continuous hour, the 7-day lifespan begins (one light for each day of therapy remaining) and continues even if unit is turned off.

See the **WARNING** in the **Alarm** section of these instructions.

Alarms

NOTE: On initial start up the green LED above the 125 mmHg selector button (default setting) or below the 75 mmHg selector button will flash until selected pressure is reached. When pressure is reached, LED will stop flashing and remain solid.

Visual Alarms - LEDs cannot be turned off or muted (paused) by the user. Visual alarms will only stop when the alarm condition has been corrected.

Audible Alarms - Two beeps which may increase in volume, depending on alarm condition. Alarms can be temporarily muted (paused) by pressing center of the alarm mute button. Beeps will re-occur after two minutes unless the alarm condition has been corrected. For additional alarm information refer to the **Alarm** section in these instructions. If alarm conditions cannot be corrected, patient should contact treating physician or caregiver.



Mute Button
(Press and Hold
center for three
seconds)

In case of an emergency contact your local emergency medical services.

WARNING: Never leave a V.A.C.® Dressing in place without active V.A.C.® Therapy for more than 2 hours. If therapy is off for more than 2 hours, remove the old dressing. Either apply a new V.A.C.® Dressing from an unopened sterile package and restart V.A.C.® Therapy; or apply an alternative dressing at the direction of the treating physician or caregiver.




Treating Physician
or Caregiver: _____ Phone: _____


For further product support see the **Contact Information** section in this guide.

Alarms

Blockage Alarm

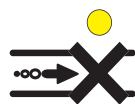

Check for Full Canister, Tubing Blockage or Closed Clamp

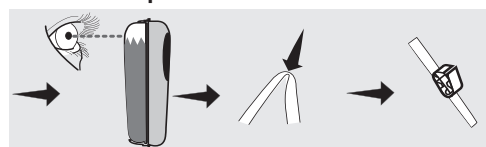


Leak Alarm




Check for Leaks


Solid yellow, therapy trying to resume.

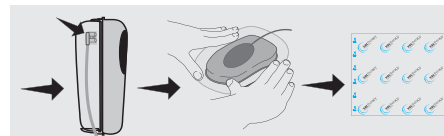


Low Battery Alarm

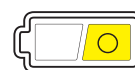


Charge Batteries

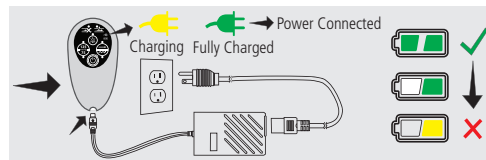

System Error Alarm

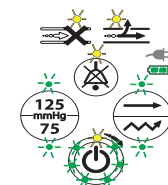

Cycle Therapy Unit OFF / ON

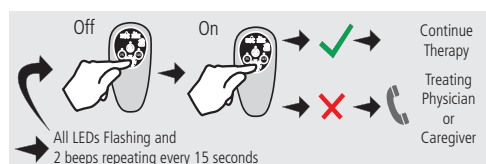

Therapy Ended Alarm



Therapy Unit will shut off



WARNING: If therapy unit shuts off, remove dressing within two hours. Replace dressing as directed by treating physician or caregiver.

Symbols Used



Conforms with the Medical Device Directive (93/42/EEC) and has been subject to the conformity procedures laid down in the council directive



ETL Listed, Conforms to UL Std. 60601-1 certified to CAN/CSA C22.2 Std. No. 601.1.



This product is designated for separate collection at an appropriate collection point. Do not dispose of as household waste.



Authorized Representative in the European Community



Consult Instructions for Use



Content Information



Catalog Number

Rx only

Federal (USA) law restricts this device to sale on or by the order of a physician.

Specifications

Environmental Conditions:

Transport and Storage Conditions:

Temperature Range: -20.2°F (-29°C) to 140°F (60°C)

Relative Humidity Range: 0-95%, non-condensing

Operating Conditions: - Temperature Range: 41°F (5°C) to 104°F (40°C)

Relative Humidity Range: 15-95%, non-condensing

Altitude Range: Optimum Performance: 0 to 8000 feet (0m to 2438m)

IP22 Protection level against ingress of solid foreign objects and liquids.

Electromagnetic Compatibility

Electromagnetic Interference - Although this equipment conforms with the intent of the directive 89/336/EEC in relation to Electromagnetic Compatibility (EMC), all electrical equipment may produce interference. If interference is suspected, move equipment away from sensitive devices or contact the manufacturer.

Radios, cell phones and similar devices may affect this equipment and should be kept at least 6.5 feet (2 meters) away from the equipment.

Manufacturer Information



Manufactured for:
KCI USA, Inc.
San Antonio, TX
78219 USA



KCI Medical Products (UK)
Ltd.
Wimborne, Dorset,
BH21 7SH
United Kingdom

Rx only



Contact Information

For additional information concerning the V.A.C. Via™ Therapy System, in the USA, call 1-800-275-4524. In other countries, contact your local KCI representative or kci-medical.com, vactherapy.com or vacvia.com.



KCI USA, Inc.
San Antonio, Texas 78219 USA
1-800-275-4524
www.kci1.com

3M™ CAVILON™ is a Trademark of 3M™ Company. All other Trademarks designated herein are proprietary to KCI Licensing, Inc., its affiliates and licensors. © 2011 KCI Licensing, Inc. All rights reserved. 360064 Rev C 09/2011

MANUFACTURER INFORMATION



Manufactured for:
KCI USA, Inc.
San Antonio, TX
78219 USA



KCI Medical Products (UK) Ltd.
Wimborne, Dorset, BH21 7SH
United Kingdom